


# CYTOGNOMIX

*Advanced reagents and software for genomic analysis*

## Chromosome Biodosimetry for Mass Radiation Emergencies

### Introduction

Ionizing radiation produces immediate, characteristic chromosome changes. Chromosome biodosimetry studies are a definitive test for occupational radiation exposure, other mass radiation emergency situations, or monitoring long term effects such as cancer. A major change results in chromosomes containing 2 constrictions, termed centromeres, instead of one (dicentric chromosomes [DCs], Figure 1). The analysis requires special expertise and microscopy equipment, is labor intensive and slow. In an emergency response to a mass casualty with exposure to a wide, unknown range of doses, biodosimetry labs will have an acute need for more methods that accurate within a few days. This proposal will streamline data acquisition using microscope system developed by Huron Technologies which will feed a software developed by Western University and Cytognomix researchers to a Huron will engineer its system to capture images 10-30 fold more quickly and Western will validate its performance with samples prepared in our lab. The project is to significantly accelerate timely radiation dose estimation, exposure or patient ascertainment, whichever occurs first.

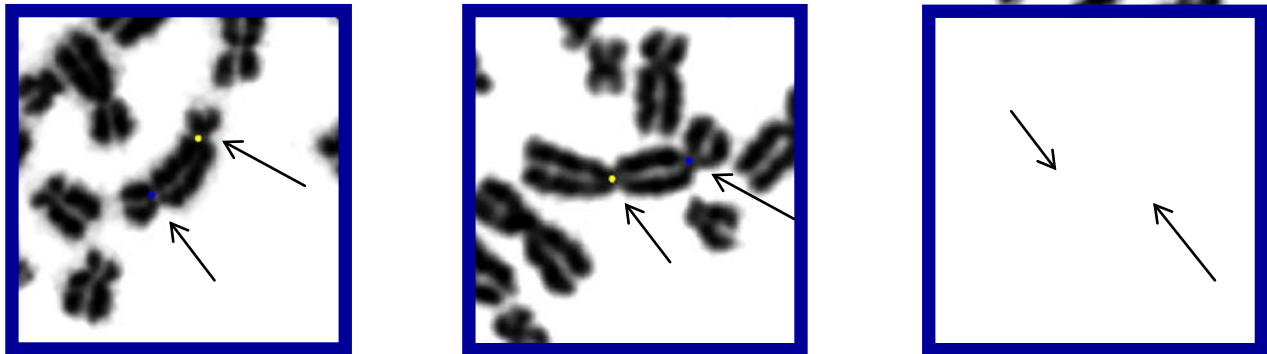


*Figure 1. Dicentric chromosomes*

### Imaging Technology

This project is a collaboration with Huron Technologies, who will carry out proof-of-principal engineering and experimental validation of a wide-field epifluorescence microscope with high power magnification capabilities ( $\geq 63X$  objective, 1.4NA). Their current wide-field TISSUEScope™ 4000 system based on the patented MACROscope technology does not have the resolution to detect dicentric chromosomes, a prerequisite for automated biodosimetry to measure radiation exposures. Huron will modify this existing system to meeting the need to generate data more quickly. Higher data throughput will maximize the benefits of the patent pending image analysis software that Western University and Cytognomix researchers are developing.

## Application Note: Chromosome Biodosimetry for Mass Radiation Emergencies

TISSUEScope™ 4000 is a third-generation fluorescence and brightfield scanner that was specifically designed for highly sensitive fluorescence imaging of up to 4 fluorophores simultaneously. The laser scan lens used in this instrument has a 20X equivalent resolution with a viewing area that is at least 10 times that of equivalent microscope objectives. This is a major advantage in whole-slide imaging, with much less potential for stitching artifacts even when imaging very large specimens. Weak and strong fluorophores can be imaged in the same scan, using adjustable gain for each channel with large dynamic range (12 bit) fluorescence detection. Figure 2 (top) shows a DAPI fluorescence image of a biodosimetry slide imaged with the TISSUEScope™ at 10X equivalent resolution. The scan area is 22mm x 46.2mm. Figure 2 (bottom) shows a 100% zoom area of the 10X image where several ROIs of interest are highlighted in green circles.

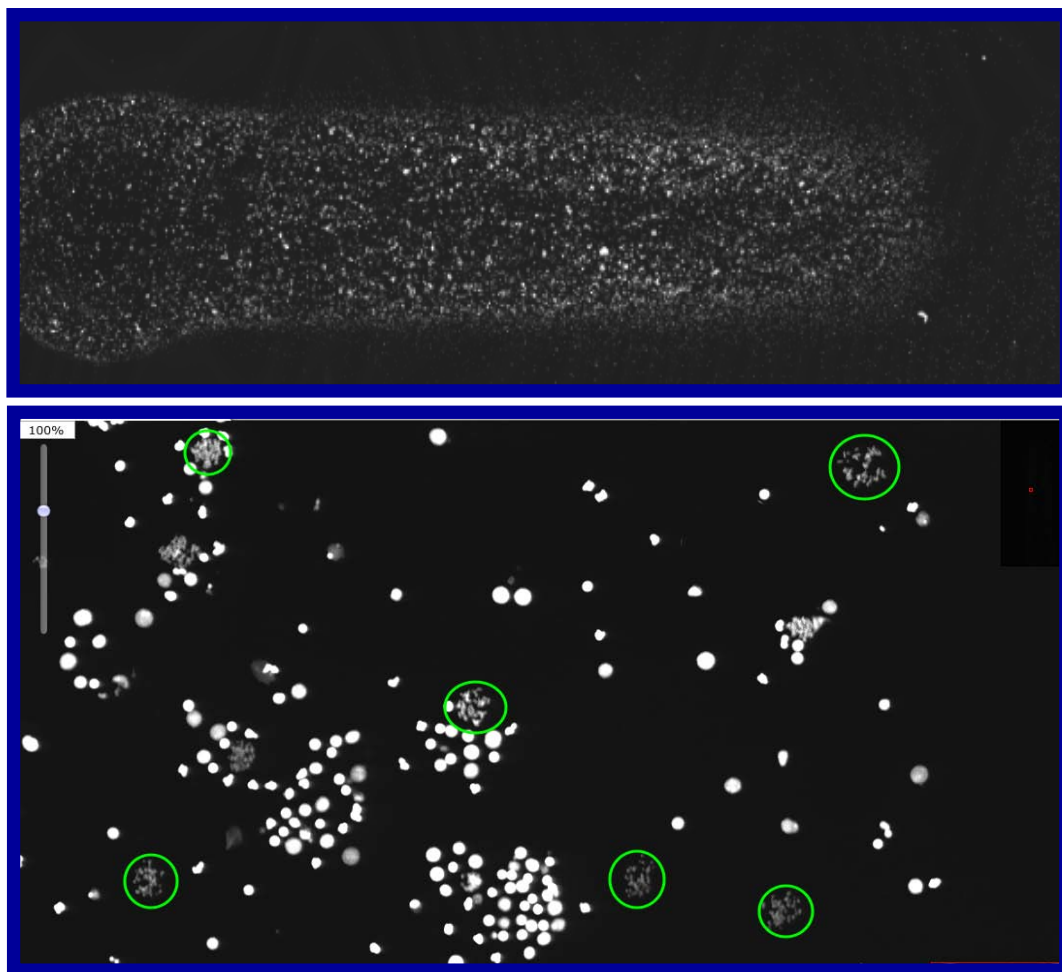


*Figure 2. (Top) 10X DAPI image of a biodosimetry slide; (Bottom) 10X image zoomed at 100% showing an area with several ROIs of metaphase cells.*

## Application Note: Chromosome Biodosimetry for Mass Radiation Emergencies

The ROIs once identified from the 10X image are imaged at 63X with a 1.4NA oil objective to clearly identify the chromosomes and using image analysis software developed by Western researchers to automatically interpret DCs. In Figure 3 images of two ROIs acquired with a microscope at 63X, at Western, showing dicentric chromosomes.

In the event of a mass casualty radiation incident, expedited preparation of samples, microscope hardware acceleration, and automated interpretation of the resulting chromosome images will be necessary to handle to overwhelming demand for biodosimetry testing. Huron's intention, in close collaboration with Cytognomix and Western University, is to develop a high speed hybrid imaging system that uses the MACROscope scanning technology and recently patented imaging technology to build a hybrid imaging system to automatically image and interpret (by implementing Cytognomix's image analysis algorithms) chromosome images. Cytognomix and Western are focusing on expediting preparation of samples and microscope hardware acceleration, and integration of these elements with automated software to interpretation of dicentric chromosome images, the latter which is currently supported by pilot funding from NIH.

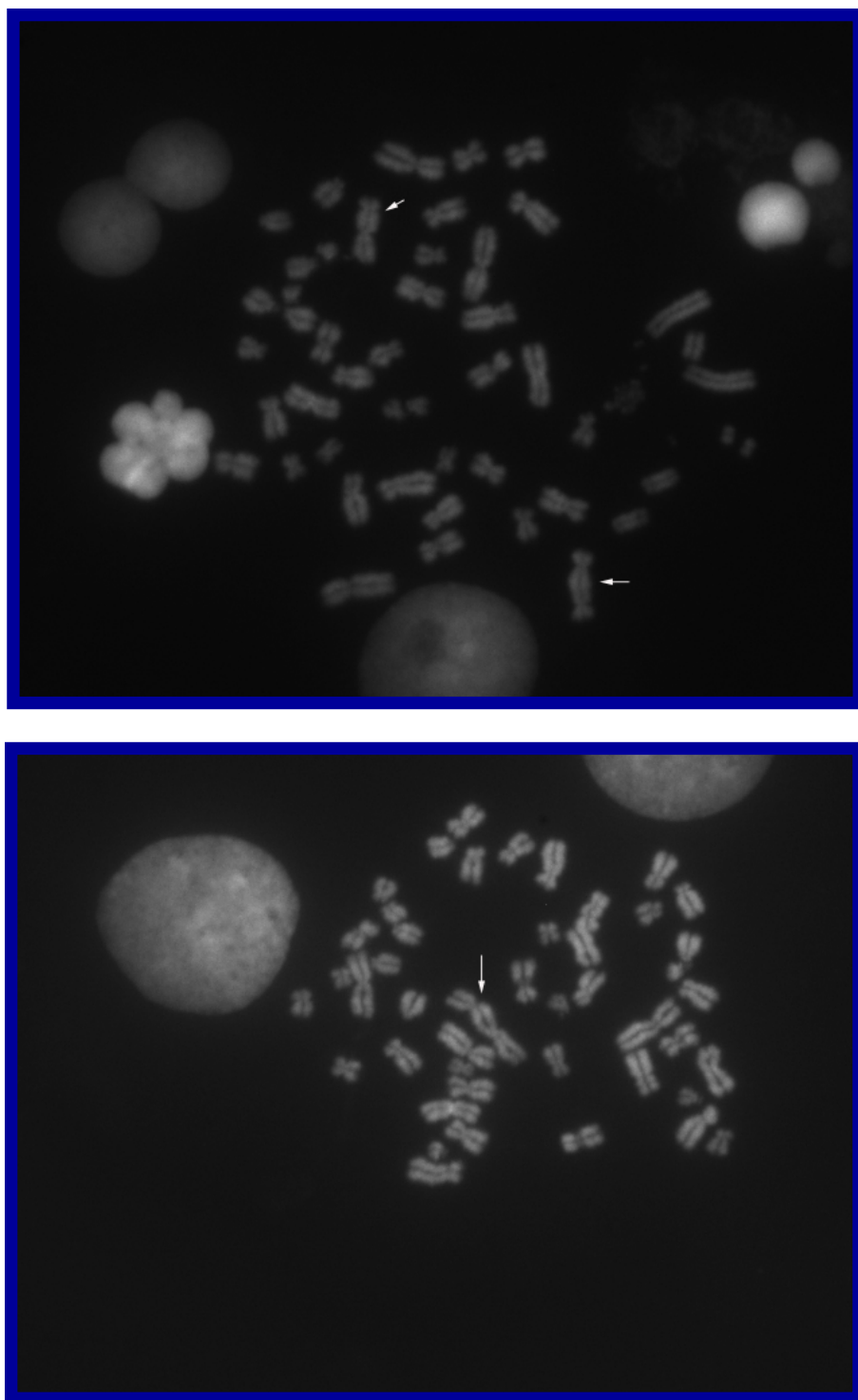

*Figure 3. Microscope images at 63X of cells with dicentric chromosomes*